


## Product Information

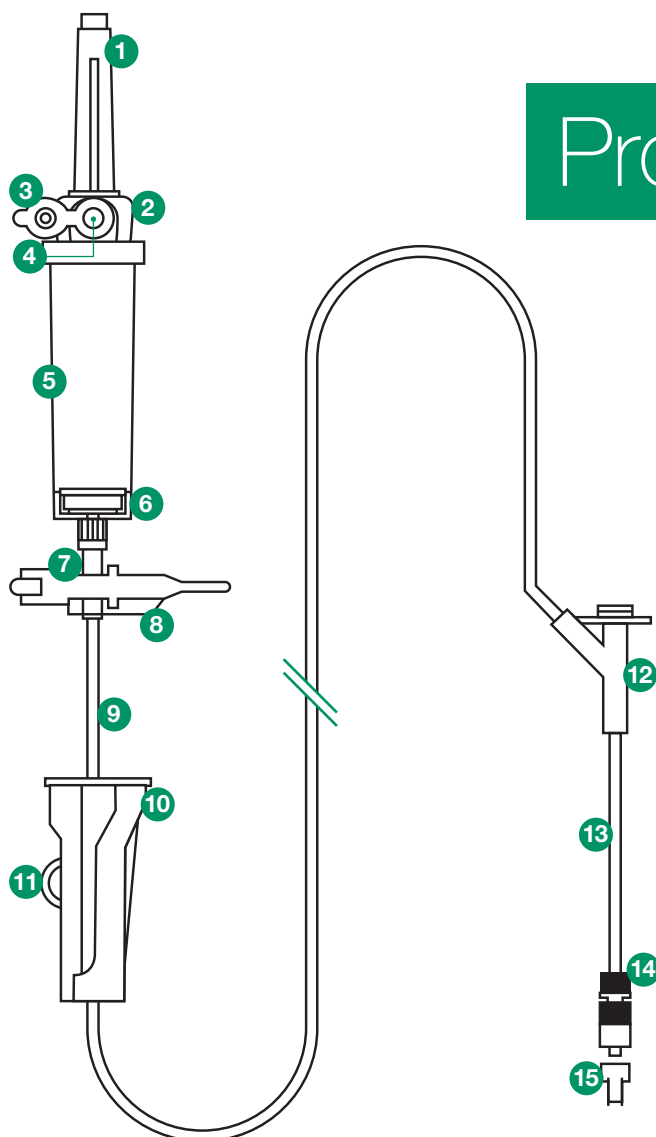

Description	Raw Materials
1 Protective cap	Polypropylene
2 Spike	ABS
3 Air filter cap	PVC
4 Air filter	Glass fiber
5 Drip chamber	PVC
6 Filter	ABS and Nylon
7 Cover tube	PVC
8 Anti-Free Flow clamp	Polyoxymethylene
9 Tubing (Ø 3.15 x 2000 mm)	PVC
10 Flow regulator housing	ABS
11 Flow regulator wheel	ABS
12 Injection point (Y-Type)	ABS + Polyisoprene
13 Tubing (Ø 3.15 x 200 mm)	PVC
14 Rotating Luer-Lock	ABS
15 Luer cap with filter	LDPE

Does not contain components made of DEHP or natural rubber latex.

## IV administration set for use with Terumo Terufusion™ infusion pumps series<sup>1</sup>



### Sterilization

With ethylene oxide to a Sterility Assurance Level of at least 10<sup>-6</sup>.



### Precautions

- Do not detach the Anti Free Flow Clip (AFF Clip ) from IV administration set.
- We recommend to switch the tubing section used in the pump after maximum 24 hr to another position as the flow rate may deviate from the set value due to deformation of the tubing section in the pump after long period of use.
- The product should be replaced and disposed of as per facility approved protocol or CDC guidelines.
- Dispose safely after use to avoid risk of infection.



## A clear view on infusion with a comprehensive portfolio of solutions.



### Packaging

#### Single Unit Pouch

#### Unit Box 25 pieces

Dimensions 275 x 230 x 135 mm, Gross Weight 1.28 kg

#### Shipping Carton 250 pieces

Dimensions 683 x 430 x 342 mm, Gross Weight 13.90 kg



### Regulatory Status

CE Certification – Class IIa medical device

**Notified Body** TÜV SÜD 0123

**Manufacturer** Poly Medicure Ltd. (India)

**EC Representative** OBELIS SA,  
Boulevard General Wahis 53, 1030 Brussels, BELGIUM

**EU Distributor** TERUMO EUROPE N.V.  
Interleuvenlaan 40, 3001 Leuven, Belgium



### Conformity

EN ISO 13485, EN ISO 11135, EN ISO 15223-1,  
EN ISO 8536-4, EN ISO 8536-8, EN ISO 10993-1,  
EN ISO 14644-1, EN ISO 11607-1, EN 868-5



### Storage

- Shelf life of 5 years
- Batch number and expiry date are printed on the label on the pouch
- Store between 5°C and 35°C



### Product Details

Spike	2-way spike
Vented	Yes
Drip chamber	20 drops/mL
Liquid filter	15 µm
Roller clamp	Orange with white wheel
Anti-Free Flow (AFF) clamp	Yes
Y-site	With needle port
Priming volume (empty drip chamber)	18.1 ml
Total length	2380 mm
Connector	Rotable lock adaptor



Does not contain components made of DEHP or natural rubber latex.

Highly accurate infusion  
combination of  
Terumo pump & sets<sup>1</sup>



**Terumo Corporation**  
+81 3 6742 8500

**Terumo Europe N.V.**  
+32 16 38 12 11

Your Terumo partner

#### Reference

<sup>1</sup> See data on file

™ Trademark  
All brand names are trademarks or registered trademarks of TERUMO CORPORATION and their respective owners.

Published by Terumo Europe N.V.